


INTRODUCING




MEDICO LEAD-LOK

Your Electrode Problem Solver

Before MEDICO LEAD-LOK

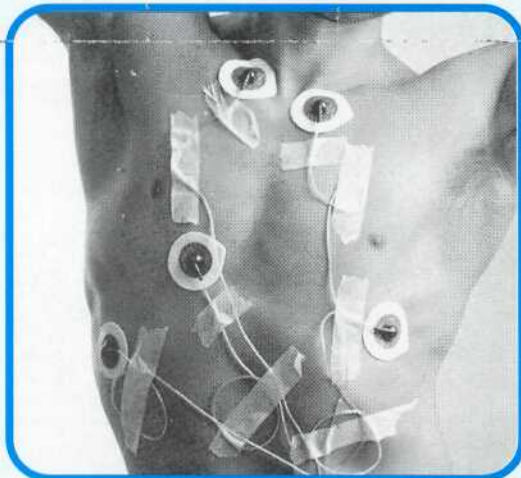

With MEDICO LEAD-LOK Model 01

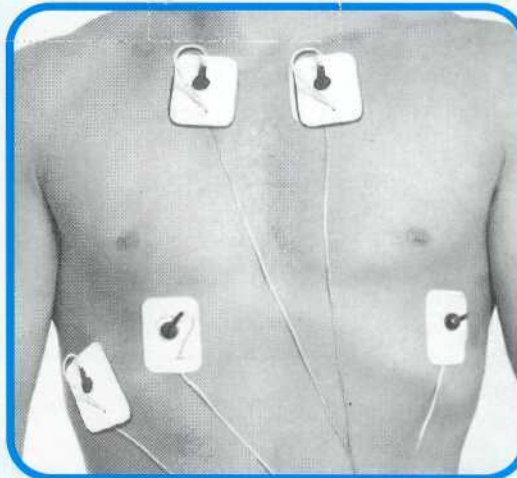

BENEFITS

- ▶ Reduces motion artifacts with lead-locking tab.
- ▶ Ensures positive, motion-free lead wire contact.
- ▶ Pregelled Liquid/Solid Gel disposable electrode.
- ▶ Biocompatible with skin tissue.
- ▶ Uses flexible, conformable high tack foam tape.
- ▶ Enables close placement when necessary.
- ▶ Rapid ECG signal recovery after Defibrillation.
- ▶ Ambulatory and long term electrode.
- ▶ Hermetically sealed packaging.
- ▶ Uses peel-off seal for impedance stability.

Backing Foam

An aggressive non-sensitizing, acrylic copolymer coated foam. Protects the skin under the electrode from liquid, air etc. thus providing for safe electrical signal passage.

Stainless Steel Stud

Strong PE Back

Label stock - free of PVC

Release Liner

Liquid or Solid Gel

Ag/AgCl Sensor

The electrical sensor is developed in such a way that the electrodes can be used together with high power reanimating equipment.

IDEAL
FOR CARDIAC
REHABILITATION

PRODUCT RANGE

MEDICO LEAD-LOK electrodes are manufactured using a special *On-Line* process. This process, in addition to being environmentally friendly, ensures that the finished product is competitively priced and is as fresh as possible.



Product code, Gel System	MLGLT-01, Liquid MSGLT-01, Solid *	MLGLT-02, Liquid	MLGLT-03, Liquid MSGLT-03, Solid	MSGLT-04, Solid	MSGLT-05, Solid
Backing, Adhesive	Foam, Aggressive	Foam, Aggressive	Foam, Aggressive	Foam, Aggressive	Foam, Aggressive
Dimensions (hxwxt)	57 x 38 x 1 mm	50 x 38 x 1 mm	50 x 48 x 1 mm	50 x 35 x 1 mm	42 x 37 x 1 mm
Suggested Use	Stress/Holter/ Echo/Monitoring	Stress/Holter/ Echo/Monitoring	Stress/Holter/ Echo/Monitoring	Stress/Holter/ Echo/Monitoring	Stress/Holter/ Echo/Monitoring
Standard Packing	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton

All Electrodes are packed 5 on a liner strip and 10 liners per pouch

* Incorporates the unique patented "Lead-Locking" tab to reduce motion artifacts



Product code, Gel System	MLGST-06, Liquid MSGST-06, Solid ‡	MLGST-07, Liquid #	MSGLT-08, Solid	MSGLT-09, Solid	MSGLT-12, Solid
Backing, Adhesive	Foam, Aggressive	Foam, Aggressive	Foam, Aggressive	Foam, Aggressive	Foam, Aggressive
Dimensions (hxwxt)	44 x 32 x 1 mm	35 x 27 x 1 mm	45 x 42 x 1 mm	53 x 34 x 1 mm	52 x 33 x 1 mm
Suggested Use	Echo/Monitoring	Echo/Monitoring	Stress/Holter/ Echo/Monitoring	Stress/Holter/ Echo/Monitoring	Stress/Holter/ Echo/Monitoring
Standard Packing	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton	50 Electrodes/pouch 10 Pouches/box 12 Boxes/carton

‡ Packed 3 Electrodes per liner strip # for Pediatric use

All Medico Lead-Lok ECG Electrodes meet or exceed the ANSI/AAMI standard EC12:2000 for disposable pregelled ECG electrodes.

MEDICO ELECTRODES (INDIA) PVT. LTD.

An ISO 9001:2000 company

t: +91 120 2424819 f: +91 120 2424859 m: +91 98104 31491

e:contact@medicoleadlok.com www.medicoleadlok.com



MEDICO LEAD-LOK



**CONTACT
US**

rev.02.so