

Control
e a s y Priming
Operation

Blood Purification Machine

Plasauto™

EZ





Features

1 *Easy Handling*

EZ «easy» *Priming* ... Priming Instructions Displayed

- Color display shows priming procedure.
- Fluid level adjustable in bubble chambers and filtrate side of plasma separator.

EZ «easy» *Control* ... All Centralized in Front

- All information is on the color display.
- All switches are located in front.

EZ «easy» *Operation* ... Touch Panel Display


- Touch panel display for set-up and data entry.

2 *For Various Applications*

- For Plasmapheresis (Plasma Exchange, Double Filtration Plasmapheresis, and Plasma Perfusion)
- Direct Hemoperfusion (DHP) is also available.


3 *Compact*

- Small size and light weight
- Retractable color display

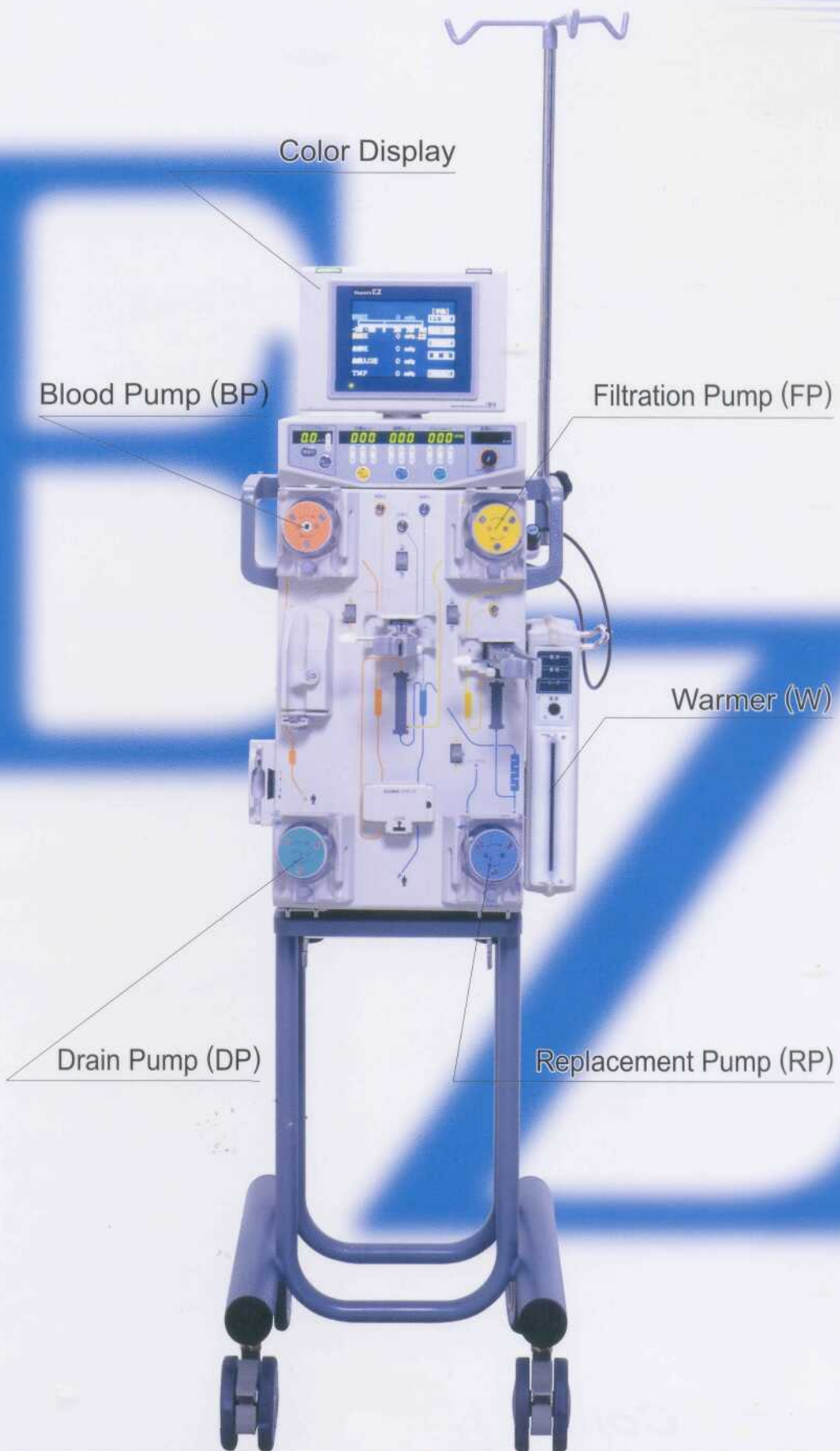


Blood Purification Machine


PE (Plasma Exchange)

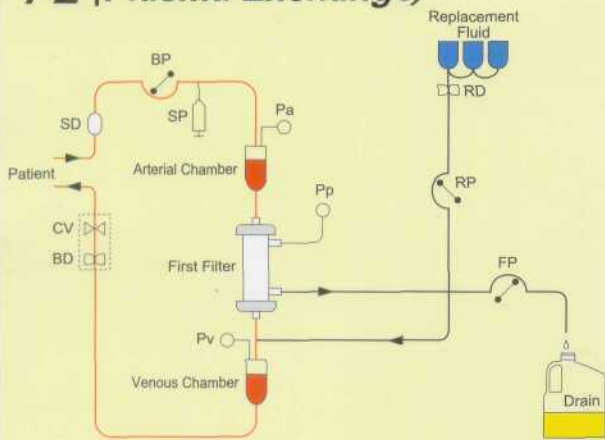

DFPP (Double Filtration Plasmapheresis)

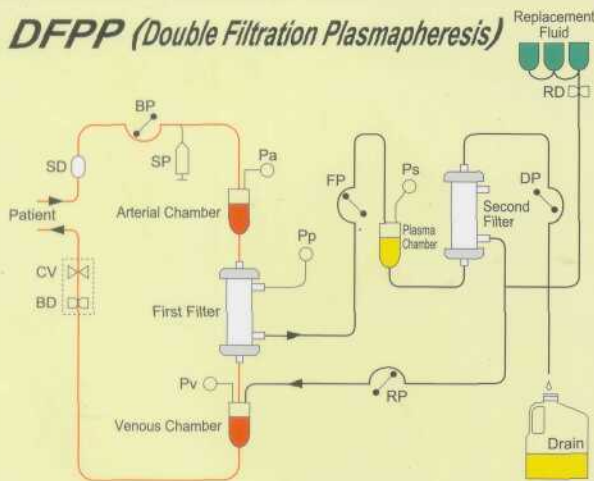

PP (Plasma Perfusion)

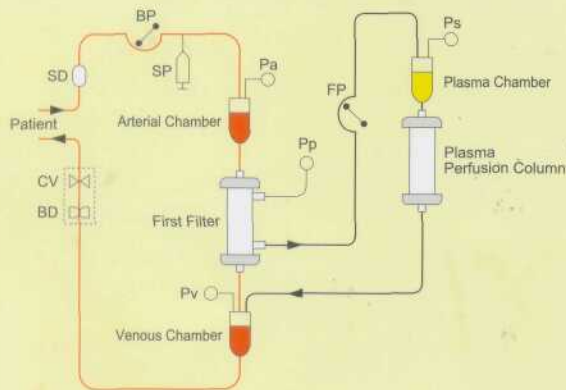

Specifications

Blood Pump (BP)

3-roller Tubing Pump
15~200mL/min.

Filtration Pump (FP), Replacement Pump (RP), Drain Pump (DP)

3-roller Tubing Pump
1~83mL/min. (0.06-5.00 L/min)

Syringe Pump (SP)

20mL Syringe
0.0~9.9mL/h

Bubble Detector (BD)

Ultrasound Type

Replacement Fluid Detector (RD)

Ultrasound Type

Clamp Valve (CV)

Electromagnetic Valve Type

Pressure Monitors

Arterial (Pa)	-400~400mmHg
Venous (Pv)	-400~400mmHg
Plasma Filtering (Pp)	-400~400mmHg
Second Filtering (Ps)	-400~400mmHg

Blood Shortage Detector (SD)

Pressure Pillow Type

Fluid Amount Monitoring

Filtration Pump Accumulation
Replacement Pump Accumulation
Drain Pump Accumulation
Blood Circulation Accumulation

Power Supply

Voltage AC110, 220V
Power Consumption 250VA

Dimension, Weight

280(W)×320(D)×1337(H)mm
Approx. 35kg



ASAHIKASEI KURARAY MEDICAL CO., LTD.

9-1, Kanda Mitoshicho, Chiyoda-ku, Tokyo 101-8482, Japan
Telephone: +81-3-3259-7724 Facsimile: +81-3-3259-7580
URL: <http://www.asahikasei-kuraray-medical.com>

Represented in Europe by:

ASAHIKASEI MEDICAL EUROPE GmbH

Arabella Center, Lyoner Straße 44-48, 60528 Frankfurt am Main, Germany
Telephone: +49-69-66371 500 Facsimile: +49-69-666 5193